BRAUN

ThermoScan®

IRT 3030

GB	English	3
CZ	Česky	13
HR	Hrvatski	23
HU	Magyar	33
PL	Polski	43
RO	Română	54
RU	Русский	65
SI	Slovenščina	77
SK	Slovenčina	87
TR	Türkçe	97
SA	اللغة العربية	116

Certain trademarks used under license from The Procter & Gamble Company or its affiliates.

ThermoScan® is a registered trademark of Helen of Troy Limited and/or its affiliates.

 Kaz Europe Sàrl
A Helen of Troy Company
Place Chauderon 18
CH-1003 Lausanne, Switzerland

© 2019, All rights reserved.

www.hot-europe.com

Designed in Germany


Patents: www.kaz.com/patents/braun

IRT3030EE

PN: A002882R0
26FEB19

Made and
printed in China

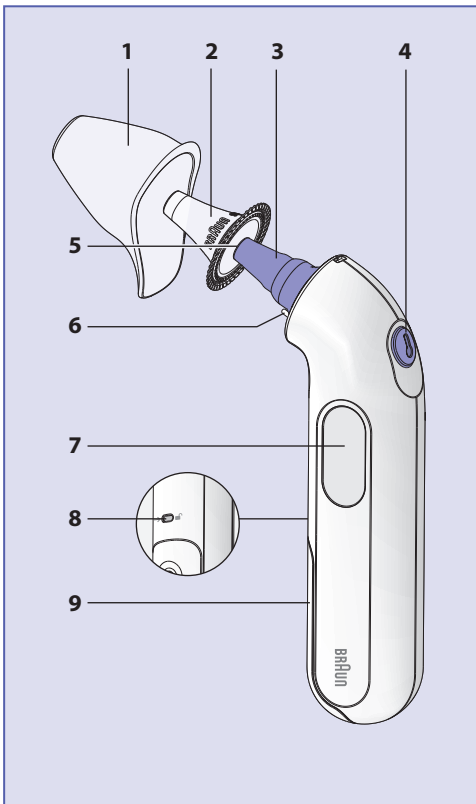Product description

- 1 Protective cap
- 2 Disposable Hygiene cap
- 3 Probe
- 4 «start» button
- 5 Probe tip (lens)
- 6 Hygiene cap detector
- 7 Display
- 8 Battery door safety lock
- 9 Battery door (1 X 3V type 2032 battery)

The Braun ThermoScan® Ear Thermometer is indicated for intermittent measurement and monitoring of human body temperature for people of all ages.

Important safety information

Please consult your doctor if you see symptoms such as unexplained irritability, vomiting, diarrhea, dehydration, changes in appetite or activity, seizure, muscle pain, shivering, stiff neck, pain when urinating, etc. in spite of absence of fever.

The operating ambient temperature range for this thermometer is 10 – 40 °C (50 – 104 °F).

Do not expose the thermometer to temperature extremes (below –25 °C / –13 °F or over 55 °C / 131 °F) or excessive humidity (> 95 % RH).

This thermometer must only be used with genuine Braun ThermoScan® disposable Hygiene caps. Other probe covers can lead to inaccuracy.

To avoid inaccurate measurements always use this thermometer with a new, clean, disposable Hygiene cap.

If the thermometer is accidentally used without a disposable Hygiene cap, clean the lens after use (see «Care and cleaning»).

Keep disposable Hygiene caps out of reach of children.

This thermometer is intended for household use only.

Use of this thermometer is not intended as a substitute for consultation with your physician. The ThermoScan® Ear thermometer is not intended to interpret hypothermic temperatures.

This thermometer must stay in stable ambient (room) temperature for 30 minutes before operating.

Audio fever indicator

Thermometer will signal by beeps normal temperature, fever or high fever.

Temperature reading	Number of beeps
Normal (≤ 37.4 °C or 99.3 °F)	One beep:)
Fever (≥ 37.5 °C or 99.5 °F)	Two beeps, repeated:))))
High fever (≥ 38.6 °C or 101.5 °F)	Four beeps, repeated:))))))))))))

Body temperature

Normal body temperature is a range. It varies by site of measurement and it tends to decrease with age. It also varies from person to person and fluctuates throughout the day. Therefore, it is important to determine normal temperature ranges.

How to use your Braun ThermoScan®

Before taking the temperature, make sure the ear is free from obstructions or excess earwax.

Do not attempt to take the reading without a disposable Hygiene cap on – the measuring sensor must be kept protected and clean from earwax to give accurate readings.

Change the disposable Hygiene cap with each use to ensure hygiene, to keep a child safe from cross-contamination and to ensure earwax residues on the disposable Hygiene cap do not impact accuracy of the reading.

- To turn the thermometer on, push the «start» button (4).


During an internal self-check, the display shows all segments.

The reading shown on the screen is the last temperature taken.


Flashing symbol of Hygiene cap will appear if no Hygiene cap is attached. Temperature can be taken only if Hygiene cap is mounted on the probe.

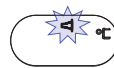
The thermometer is ready to use after two beeps are sounded and the non-flashing ear icon appears.

- In clinical studies, precision and accuracy of ear thermometer measurement has been shown to be influenced by variability in user technique. Correct probe placement is essential for accurate measurements. Place the ear probe snugly in ear canal and direct towards opposite temple, then push the "start" button (4).

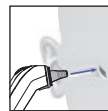Keep the thermometer steady in the ear canal.

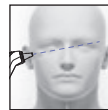Temperature in right and left ear differs, always compare measurements from the same ear. The natural physiological difference may reach 0.2 °C. In the absence of confounding factors, in high-risk population (like infants or immunodepressed people) use the highest of the two to interpret the patient's condition.

One beep signals normal temperature. Audio fever indicator will trigger 2 consecutive beeps when the reading is equal or higher to 37.5 °C (99.5 °F) and 4 beeps if the reading is equal or higher to 38.6 °C (101.5 °F). The result is shown on the display (7).

3. Gently pull off the disposable Hygiene cap (2).

Do not use force in removing the disposable Hygiene cap.


4. Put on a new, clean, disposable Hygiene cap by snapping into position.

5. Place the protective cap on (1).

The Braun ThermoScan Ear thermometer turns off automatically after 60 seconds of inactivity. The thermometer can be turned off by pressing and holding the start button for more than 5 seconds.

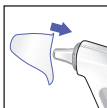Disposable Hygiene caps

Use a clean disposable Hygiene cap with each measurement:

1. To ensure hygienic use
2. To protect your child from cross-contamination
3. To ensure earwax residues do not impact accuracy of the reading.

Earwax and dirt on the tip of the disposable Hygiene cap, even in microscopic layers, can create an obstacle between the sensor and the infrared heat emitted by the ear and impact the accuracy of the measurement.

Temperature-taking hints


- Replace the disposable Hygiene cap with every use to maintain accuracy, hygiene and to prevent cross-contamination.
- A measurement taken in the right ear may differ from the measurement taken in the left ear. Therefore, always take the temperature in the same ear.
- The ear must be free from obstructions or excess earwax build-up to take an accurate reading.
- External factors may influence ear temperatures, including when an individual has:
 - been lying on one ear or the other
 - had their ears covered
 - been exposed to very hot or very cold temperatures, or
 - been recently swimming or bathing.
 In these cases, remove the individual from the situation and wait 20 minutes prior to taking a temperature.
- If ear drops or other ear medications have been placed in the ear canal of one ear, take the temperature in the untreated ear.

Care and cleaning

The probe tip (5) is the most delicate part of the thermometer. It has to be clean and intact to ensure accurate readings.

If the thermometer is accidentally used without a disposable Hygiene cap, clean the probe tip as follows:

Very gently wipe the surface with a cotton swab or soft cloth moistened with alcohol. Do not use cleaners other than alcohol on the probe tip. Use of other cleaners could damage the tip.

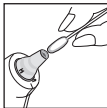Once the alcohol has dried completely, a new disposable Hygiene cap can be put on and a temperature measurement taken. If the probe tip is damaged, please contact your authorized local service center.

Use a soft, dry cloth to clean the thermometer display and product exterior. Do not use abrasive cleaners. Never submerge this thermometer in water or any other liquid.

Store the thermometer and disposable Hygiene caps in a dry location free from dust and contamination and away from direct sunlight.

Additional disposable Hygiene caps (LF40) are available at most stores selling Braun ThermoScan.

Trouble-shooting

Situation	Solution
 No disposable Hygiene cap is attached.	Attach new, clean, disposable Hygiene cap.
 Measurement before device stabilization.	Wait until all the icons stop flashing.
 Ambient temperature is not within the allowed operating range (10–40 °C or 50–104 °F).	Allow the thermometer to remain for 30 minutes in a room where the temperature is between 10 and 40 °C or 50 and 104 °F.

	Temperature taken is not within typical human temperature range (34–42.2 °C or 93.2–108 °F).	Make sure the probe tip and lens are clean and a new, clean, disposable Hygiene cap is attached. Make sure the thermometer is properly inserted. Then, take a new temperature.
		
	System error.	Reset the thermometer by removing the battery and then replacing it again, wait 1 minute, then turn on again. Please contact your authorized Service Center (See Guarantee Card).
	Battery is low.	Insert new battery.
	Battery is too low. Device cannot be powered on to the ready state.	Insert new battery.
Do you have any further questions?		Please contact your authorized Service Center (See Guarantee Card).

Changing the temperature scale

Make sure the thermometer is turned off. Press and hold down the <<start>> button. Keep holding down the <<start>> button and after about 3 seconds the display will start to cycle between temperature scales (°C or °F).


Release the start button when the desired scale is shown. This will reset the thermometer to the new temperature scale.

Replacing the battery

The thermometer is supplied with one 3V Lithium cell battery (type 2032). Insert new battery when the low battery symbol appears on the display.

The battery door design is intended to prevent a small child from easily opening the battery door to address a possible choke hazard. Open the battery compartment by inserting a narrow pointed object such as a ball point pen or paperclip into the small hole above the battery door (8). While pressing into the hole slowly slide back the battery door with the other hand.

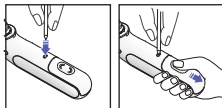To close battery door align the tabs along the side of the battery into the battery compartment and slowly slide back into position until the door clicks into its locked position.

Calibration

This product has been calibrated at the time of manufacture. If used according to the use instructions, periodic calibration is not required. This device is not intended as a substitute for regular check-ups by your doctor, please consult your doctor if you have any doubt about the temperature reading. LOT and Serial numbers of your product must be provided as they are essential to record and follow up your inquiry or claim. The LOT and Serial numbers are located in the battery compartment.

The production date is given by the LOT number and can be deciphered as explained below:

The first 3 numeric digits in the LOT number represent the day of the year of manufacture. The next 2 numeric digits represent the last two numbers of the calendar year of manufacture and the letter(s) at the end designate the manufacturer of the product. (E.g.: LOT No.: 12313tav this product was made on the day 123, year 2013 at the manufacturer code tav.)

Guarantee

EMC tables and Consumer Card are available on our website at www.hot-europe.com/support.

For support requests, please visit www.hot-europe.com/support or find service contact information on the last page of this manual.

UK only: This guarantee in no way affects your rights under statutory law.

Product specifications

Model:	IRT 3030
Displayed temperature range:	34 °C – 42.2 °C (93.2 °F – 108 °F)
Operating ambient temperature range:	10 °C – 40 °C (50 °F – 104 °F)
Storage temperature:	-25 °C – 55 °C
Display resolution:	0.1 °C or °F
Accuracy for displayed temperature range:	± 0.2 °C (35.5 – 42 °C) (95.9 – 107.6 °F) ± 0.3 °C (outside this temperature range)
Clinical bias:	-0.11 °C (-0.198 °F)
Limits of agreement:	± 0.25 °C (0.45 °F)
Clinical repeatability:	± 0.10 °C (0.18 °F)
Measuring Site:	Ear
Reference Site:	Oral (estimated oral temperature is displayed)
Operating and storage relative humidity:	≤ 95% non-condensing
Battery life:	2 years /1000 measurements
Service life:	5 years
Warranty:	2 years

If device is not used or stored within specified temperature and humidity ranges the technical accuracy cannot be ensured.

Type BF
equipment

See instructions
for use

Operating temperature

Storage temperature

Subject to change without notice.

IP22: Protected against solid foreign objects of 12.5 mm diameter and greater. Protected against vertically falling water drops when the device is tilted up to 15 degrees.

EN 12470-5:2003 Clinical thermometers—Part 5: Performance of infra-red ear thermometers (with maximum device).

ISO 80601-2-56:2009 Medical electrical equipment—Part 2-56: Particular requirements for basic safety and essential performance of clinical thermometers for body temperature measurement.

This product conforms to the provisions of the EC directives 93/42/EEC (Medical Device Directive) and 2011/65/EU (RoHS).

The Braun ThermoScan needs special precautions regarding EMC. For detailed description of EMC requirements please refer to www.hot-europe.com/support.

Portable and mobile RF communications equipment can affect the Braun ThermoScan.

This product contains batteries and recyclable electronic waste. To protect the environment, do not dispose of it in the household waste, but take it to appropriate local collection points.

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