tria®

HAIR REMOVAL LASER 4X

Quick Start Guide

USEFUL INFORMATION

- Hair follicles that have been deactivated by the Tria Laser 4X will not grow back.
- Laser hair removal does not cause hair to fall out immediately. After your first treatment, there will be no visible change to the hair.
- Some users may feel the sensation of a rubber band snapping on the skin. Others will feel a warm, prickling sensation. As the amount of hair decreases, the sensation will also decrease.
- In a clinical study, higher levels were found to be 2-3 times more effective. So for best results, we recommend you use the highest setting that you can tolerate.
- Treat every 2 weeks for up to 3 months or until you reach the desired results.

Clinical studies on file

CHECK YOUR SKIN TONE

The Tria Laser 4X is intended only for light to medium skin tones and medium to dark hair. It does not work on white, grey, blonde or red hair. If your skin tone matches any of the skin tones labeled "Not Safe" for treatment, then DO NOT USE the Tria Laser 4X.

SKIN TONE CHART

INTENDED SKIN TONES:

NOT SAFE:

Tria Beauty Customer Care
United States
877-321-8742
customer@triabeauty.com
triabeauty.com


tria®

 Tria Beauty, a company of Channel Investments, LLC
7999 Gateway Boulevard, Suite 100, Newark, CA 94560 USA
All rights reserved. Tria, Tria Beauty, and the Tria logo are trademarks of Tria Beauty.
©2017

15017B

GET TO KNOW YOUR TRIA LASER 4X

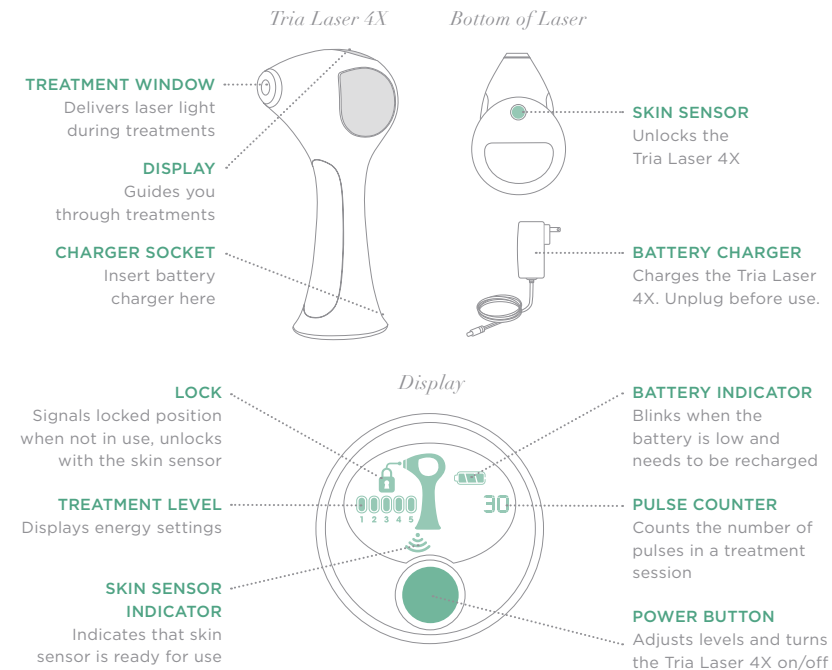BEFORE FIRST TREATMENT

For detailed directions and important safety warnings, please read the complete *Instructions for Use* guide before use.

CHARGE: Charge the Tria Laser 4X until the battery light is full (~2 hours). Unplug before use.

ACTIVATE: A one-time, 5-minute activation is required before your first treatment. In the US: visit triabeauty.com/activate or call Customer Care at 877-321-8742. In Canada: visit triabeauty.ca/activate or call Customer Services at 877-321-8742.

HOW TO TREAT

Before every treatment, shave, clean and dry the area you wish to treat. Before your first full treatment on any new area, perform a patch test in that area. Repeat treatments every 2 weeks for up to 3 months or until you're satisfied with results.


1 TURN ON DEVICE

Unplug the Tria Laser 4X from the battery charger and press the power button.


2 UNLOCK

Unlock the device by holding the Skin Sensor against the skin in the area you wish to treat. If the Tria Laser 4X is safe for your skin tone, a tone will sound and the Tria Laser 4X will unlock (if device doesn't unlock, see *Instructions for Use* for more information).

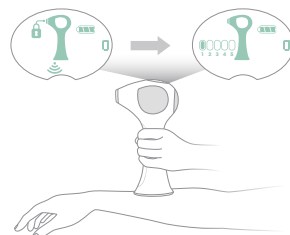
3 SET THE TREATMENT LEVEL

Press the power button to adjust levels. For best results, use the highest setting that you can tolerate.

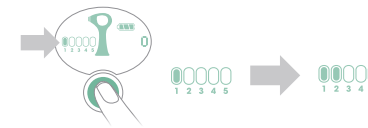
4 TREAT. HOLD FOR 2 BEEPS

- Touch and hold the treatment window to the skin until you hear two beeps.
- After the two beeps have sounded, lift the laser and place it back on the skin so that it overlaps the previous spot (about 1/4 inch or 5 mm).

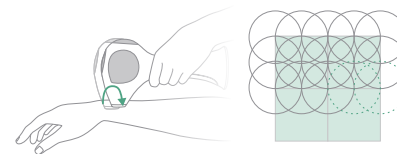5 KEEP GOING TO COVER THE AREA

Use the Tria Laser 4X Pulse Counter as a guide to how many pulses you may need to cover each area per treatment.

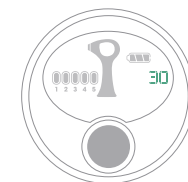PULSE COUNTER GUIDE

UPPER LIP
Pulses: 25
Minutes: 2-3

UNDERARMS
Pulses: 100 per side
Minutes: 4-5 per side

BIKINI
Pulses: 200 per side
Minutes: 8-10 per side

UPPER LEG
Pulses: 600 per side
Minutes: 25-30 per side

LOWER LEG
Pulses: 600 per side
Minutes: 25-30 per side

Pulse counts and times are only approximations—your actual treatment area and number of pulses may vary. Always follow the treatment instructions and warnings provided in the *Instructions for Use* booklet. When treating larger areas like legs and back, you may need more than 600 pulses to cover the area.

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