HAIR REMOVAL LASER 4X Quick Start Guide

CHECK YOUR SKIN TONE

USEFUL INFORMATION

Hair follicles that have been deactivated

by the Tria Laser 4X will not grow back.

 Laser hair removal does not cause hair to fall out immediately. After your first treatment, there will be no visible change

Some users may feel the sensation of a

rubber band snapping on the skin. Others

amount of hair decreases, the sensation

• In a clinical study, higher levels were found

• Treat every 2 weeks for up to 3 months or

until you reach the desired results.

setting that you can tolerate.

to be 2-3 times more effective. So for best results, we recommend you use the highest

will feel a warm, prickling sensation. As the

to the hair.

will also decrease.

The Tria Laser 4X is intended only for light to medium skin tones and medium to dark hair. It does not work on white, grey, blonde or red hair. If your skin tone matches any of the skin tones labeled "Not Safe" for treatment, then DO NOT USE the Tria Laser 4X.

SKIN TONE CHART

INTENDED SKIN TONES:

Brown or Black

NOT SAFE:

Tria Beauty Customer Care 877-321-8742

customercare@triabeauty.com triabeauty.com


Tria Beauty, a company of Channel Investments, LLC 7999 Gateway Boulevard, Suite 100, Newark, CA 94560 USA All rights reserved. Tria, Tria Beauty, and the Tria logo are trademarks of Tria Beauty.

15017B

Clinical studies on file

GET TO KNOW YOUR TRIA LASER 4X

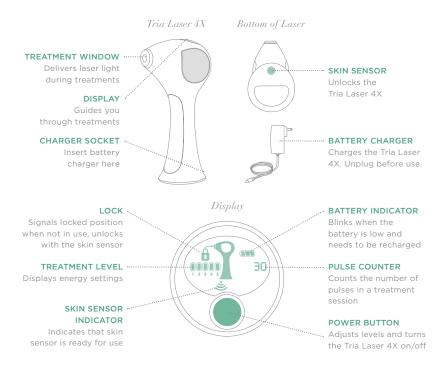BEFORE FIRST TREATMEN

For detailed directions and important safety warnings, please read the complete *Instructions for Use* guide before use.

HOW TO TREAT

Before every treatment, shave, clean and dry the area you wish to treat.

Before your first full treatment on any new area, perform a patch test in that area.

Repeat treatments every 2 weeks for up to 3 months or until you're satisfied with results.


TURN ON DEVICE

Unplug the Tria Laser 4X from the battery charger and press the power button.

2 UNLOCK

Unlock the device by holding the Skin Sensor against the skin in the area you wish to treat. If the Tria Laser 4X is safe for your skin tone, a tone will sound and the Tria Laser 4X will unlock (if device doesn't unlock, see *Instructions for Use* for more information).

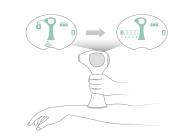CHARGE: Charge the Tria Laser 4X until the battery light is full (-2 hours). Unplug before use

CTIVATE: A one-time, 5-minute activation is required before your first treatment. the US: visit triabeauty.com/activate or call Customer Care at 877-321-8742. Canada: visit triabeauty.ca/activate or call Customer Services at 877-321-8742.

3 SET THE TREATMENT LEVEL


Press the power button to adjust levels.

For best results, use the highest setting that you can tolerate.

TREAT. HOLD FOR 2 BEEPS

- Touch and hold the treatment window to the skin until you hear two beeps.
- After the two beeps have sounded, lift the laser and place it back on the skin so that it overlaps the previous spot (about 1/4 inch or 5 mm).

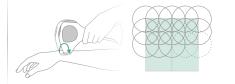
5 KEEP GOING TO COVER THE AREA

Use the Tria Laser 4X Pulse Counter as a guide to how many pulses you may need to cover each area per treatment.

PULSE COUNTER GUIDE

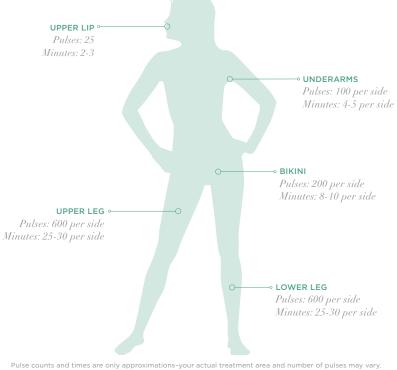ulse counts and times are only approximations-your actual treatment area and number of pulses may var Always follow the treatment instructions and warnings provided in the Instructions for Use booklet. When treating larger areas like legs and back, you may need more than 600 pulses to cover the area.

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